

# HümiOne™

## DISPOSABLE PREFILLED BUBBLE HUMIDIFIER



Single Use

High Quality Disposable  
Prefilled Bubble Humidifier

USP Compliance  
Sterile Water

Individual Packing  
With Adaptor

Clinically Proven  
Simple Yet Practical



STERILE

SBS

MD



CE 2460

UDI

Imported by :

**Oneness Enterprises**  
Medical Equipments / Surgical Disposables  
Sales & Services

Marketed by :

**PureCure+**  
Biomedical Devices Pvt. Ltd.

Shop No. 1, Galaxy Corner, Near Kailas Jeevan Factory,  
Sai Puram Busstop, Dhayari, Pune 41.

+91 8600005800 / 9766939741

onenessenterprises15@gmail.com

## Product Intended Use:

The intended use of this device is to moisturize breathing gases for patients. HumiOne is designed to be used with various breathing gas sources including oxygen concentrators, gas cylinders and wall outlets which are connected to patients via cannulas and face masks. HumiOne is also highly recommended for hazardous patients or isolation wards to help eliminate the risk of cross-contamination.



### CDH350

- Volume 350ml
- Weight 390g
- Flow Rate 0-15 LPM



### CDH340

- Volume 340ml
- Weight 380g
- Flow Rate 0-15 LPM



### CDH450

- Volume 450ml
- Weight 490g
- Flow Rate 0-6 LPM



### CDH500s

- Volume 500ml
- Weight 540g
- Flow Rate 0-6 LPM

Oxygen Pressure ( Range ) As low as 0.5 bar ( @ 5 LPM )

## SAFETY PRECAUTION

- Single use only
- Non - pyrogenic sterilized water.
- DO NOT USE if seal is tempered.
- Discard used bottle, Do not reuse.

## USAGE

- Use immediately after nozzle tab is pulled and removed.

## OTHER INFORMATION

- Universal Male Adaptor for Oxygen Tubing at Nozzle .
- Universal Adaptor for Oxygen Flow Meter at Connector .
- Side bottle volume scales for easy monitoring .

## SPECIFICATIONS


### Material

Bottle : Polypropylene (PP)  
Adaptor : Polystyrene (PS) and Polypropylene (PP)


Water : Sterile water for inhalation. Complies to United State Pharmacopoeia (USP) for bodily safe, non-pyrogenic and effective humidity for inhalation of dry medical oxygen

Dimension : W 95mm x D 85mm x H 150mm  
Box Dimension : 55cm (L) x 43cm (W) x 20cm (H)  
Expiry Date : 5 Years after manufacturing date

## INSTRUCTIONS FOR USE

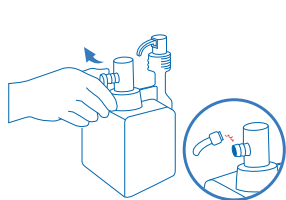


1. PULL the nozzle tab upwards and discard.  
**DO NOT TWIST!**

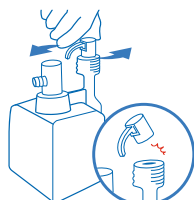


2. TWIST the air flow tab and discard.




3. Fasten the connector to HumiOne™

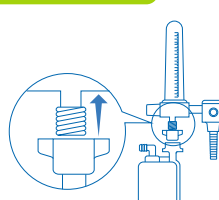

4. Mount connector to the oxygen outlet.

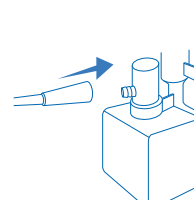

5. Connect the oxygen tubing.

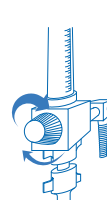

6. Turn on the oxygen supply.

**HumiOne™**

Should form bubbles. Check for leakages.